



VISCOELASTIC MEDICAL DEVICES




Siloam Ophthalmics

No:23, Level 1, Periyar Street, Kannikapuram,
Avadi, Chennai – 600054 , Tamil Nadu, INDIA.

✉ info@rainbowoculardevices.com ☎ +60 10-266 1852

VEU-VISC 2ML PFS

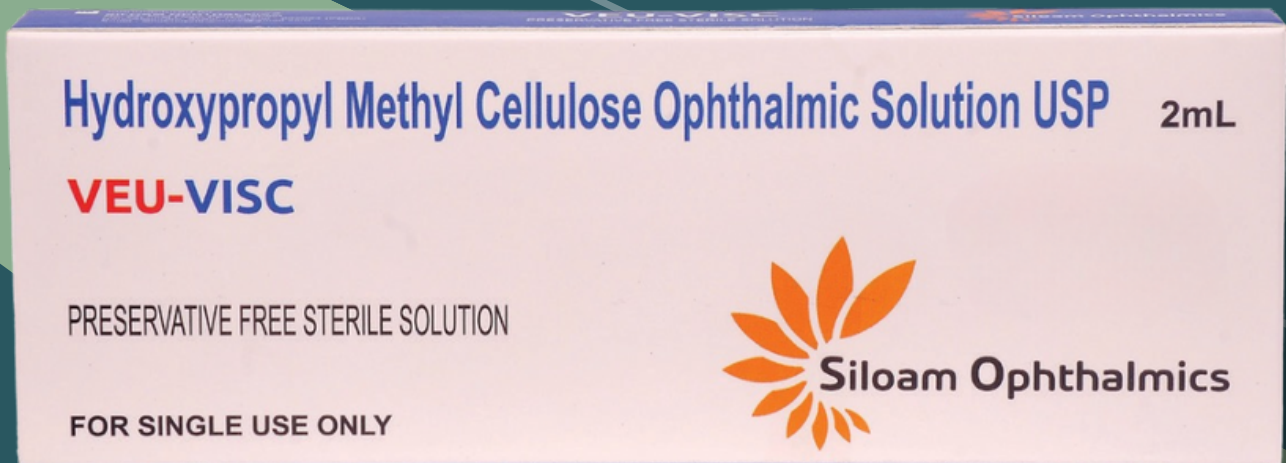

PRODUCT DESCRIPTION

Veu-Visc PFS is a viscoelastic solution of optimum viscosity. It contains highly purified grade of HPMC 2%w/v. It is clear isotonic sterile, non-inflammatory, used as a lubricant and protector in IOL surgery.

Composition

Each mL contains:
Hydroxy Propyl Methyl
Cellulose USP ...2%w/v
Sterile Isotonic base q.s.

Packaging:

Primary Packing:- 3 ml
syringes, plunger stopper,
plunger rod & pouch
Secondary Packing:- Carton
box with 23G Sterile Cannula .

Shelf Life

24 months

Storage:

Store between 15°C-25°C

Product Features:

- Clear Visualization
- Preservative free sterile solution
- Optimum Viscosity
- Better chamber maintenance
- Friendly with corneal endothelium
- Provide excellent Lubrication
- Simple to remove from the anterior chamber

Dosage & Administration

VEU-VISC (Hydroxypropyl Methyl Cellulose Ophthalmic Solution USP) should be carefully introduced into the anterior chamber using the sterile cannula. It may be applied into the chamber prior to or following delivery of intraocular lens (IOL). It may be used to coat an intraocular lens as well as tips of surgical instruments prior to implantation which protect the corneal endothelium from possible damage during cataract surgical procedure.


Indication

VEU-VISC is indicated as a surgical aid (medical device) in anterior segment surgical procedures involving the anterior chamber of the eye, including extraction of the lens and insertion of intraocular lenses. It maintains the depth of the anterior chamber during the whole surgical procedure and permits greater operative precision without the risk of damaging the endothelium of the cornea or other intraocular tissues.

Supply

VEU-VISC (Hydroxy Propyl Methyl Cellulose Ophthalmic Solution USP 2%w/v) is a sterile, non- pyrogenic, viscoelastic preparation supplied in 3mL syringe with luer tip and luer lock packed in double pouch and terminally sterilized.

VEU-VISC PLUS 1.4

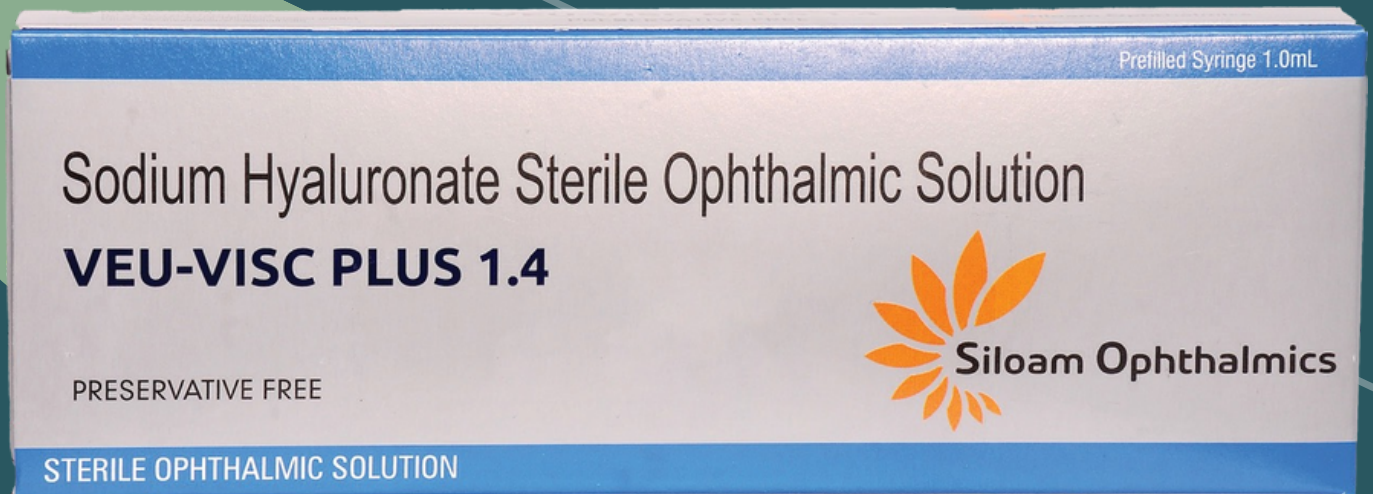

PRODUCT DESCRIPTION

VUE-VISC PLUS 1.4 ophthalmic visco surgical device is a sterile, non-pyrogenic, high molecular weight, Non- inflammatory highly purified fraction of Sodium Hyaluronate dissolved in physiological Sodium Chloride, Phosphate buffer. Each mL of Veu-Visc Plus contains 14mg of Sodium Hyaluronate in a sterile isotonic solution.

Composition

Each mL contains:
Sodium Hyaluronate BP
Sterile Isotonic Solution
14mg q.s.

Packaging:

Primary pack:– 2.25ml glass
Syringe & pouch.
Secondary package :–Carton
box with 27G Sterile Cannula.

Shelf Life

24 months

Storage:

Store between 8°C-25°C



Direction of use

VUE-VISC PLUS 1.4 Cataract surgery and intraocular lens (IOL) implantation: Instill VUE-VISC PLUS 1.4 slowly in to the artificially created opening in the eye through the blunt cannula into the anterior chamber. The use of VUE-VISC PLUS 1.4 is most effective when the Instill is made before phacoemulsification, removal of the cataract and before insertion of the IOL. VUE-VISC PLUS 1.4 can also be applied to the IOL before placement. During the procedure, more VUE-VISC PLUS 1.4 can be infused for anterior chamber maintenance or to replace viscoelastic lost during surgery.

Indication


VUE-VISC PLUS 1.4 is viscoelastic ophthalmic devices indicated as an aid for ophthalmic surgical procedures and anterior segment during cataract extraction and IOL implantation. VUE-VISC PLUS 1.4 helps to maintain deep anterior chamber during anterior segment surgery reduce trauma to the corneal endothelium and surrounding ocular tissue It helps to push back the vitreous face and prevent formation of a flat chamber during the surgery.

Supply

VUE-VISC PLUS 1.4 is a sterile, single use viscoelastic device supplied in a non-pyrogenic Glass Syringe containing 1ml & 2mL with single use sterile angular cannula 27 G.



VEU-VISC PLUS 1.8

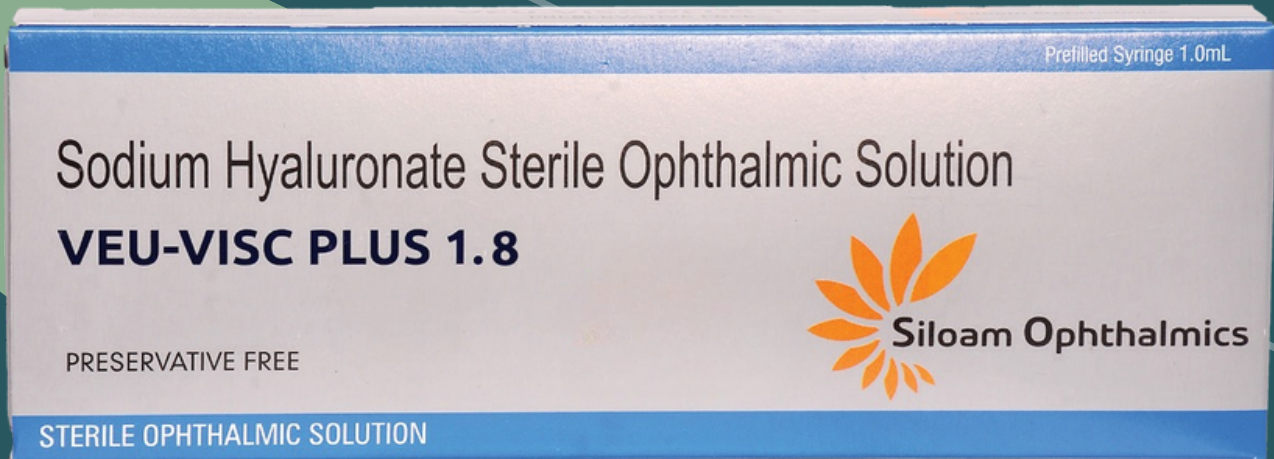

PRODUCT DESCRIPTION

VEU-VISC Plus 1.8 ophthalmic visco surgical device is a sterile, non-pyrogenic, high molecular weight, Non-inflammatory highly purified fraction of Sodium Hyaluronate dissolved in physiological Sodium Chloride, Phosphate buffer. Each mL of VUE-VISC contains 18mg of Sodium Hyaluronate in a sterile isotonic solution.

Composition

Each mL contains:
Sodium Hyaluronate BP
Sterile Isotonic Solution
18mg q.s.

Packaging:

Primary pack :- 2.25mL glass
Syringe & pouch Secondary
package :- Carton box with
27G Sterile Cannula`

Shelf Life

24 months

Storage:

Store between 8°C-25°C



CONTACT US



Siloam Ophthalmics

**No:23, Level 1, Periyar Street,
Kannikapuram, Avadi, Chennai – 600054 ,
Tamil Nadu, INDIA.**

Malaysia Distributor

**78 A, Jalan BP 7/8, Bandar Bukit Puchong,
47120 Puchong, Selangor, MALAYSIA**

 +60-3-8066-7448

 +60 10-266 1852

+60 10-240 5768

 www.rainbowoculardevices.com

 info@rainbowoculardevices.com